

Controller | QMS 800



QC 800 – QUADROPLE CONTROLLER

New OPC-UA communication protocol designed into latest generation.



IS 816 – ION SUPPLY

The ion source is controlled by the IS 816, with all parameters defined and saved in the module.



Dimensions: 400 mm L x 483 mm W x 133 mm H



IO 821 – I/O SYSTEM

Allows connection of analog and digital inputs and outputs, such as an external total pressure gauge.





HV 801 – HIGH-VOLTAGE SUPPLY

Supplies the voltage required to operate the SEM 217 and supports voltage signal communication between the QMS and QMH.

QMG 800

| SPECIFICATIONS | QMA 430 | QMA 400 |
|--|-------------------------------------|-------------|
| Mass Range | 1-300 AMU | 1-512 AMU |
| Rod Material | Stainless Steel | Molybdenum |
| Maximum Operating Pressure | 1E-5 mbar | |
| Faraday Sensitivity | | |
| Axial Ion Source | 2E-4 A/mbar | 3E-4 A/mbar |
| Cross Beam Ion Source | 8E-5 A/mbar | 7E-5 A/mbar |
| Cross Beam Gas Tight Ion Source | 8E-5 A/mbar | 7E-5 A/mbar |
| Cross Beam Ion Source with Magnets | 7E-4 A/mbar | 7E-4 A/mbar |
| Cross Beam Gas Tight Ion Source with Magnets | 7E-4 A/mbar | 7E-4 A/mbar |
| Grid Ion Source | 2E-4 A/mbar | 2E-4 A/mbar |
| Resolution | 0.3 – 7.0 AMU | |
| Minimum Detectable Partial Pressure | 2E-15 mbar | 1E-15 mbar |
| Partial Pressure Ratio with Electron Multiplier | < 1 ppb | < 0.5 ppb |
| Maximum Operating Temperature (Sensor & Electronics) | Sensor: 150°C Electronics: 400°C | |
| Measurement Speed | 125 µs | |

CONFIGURATION OPTIONS

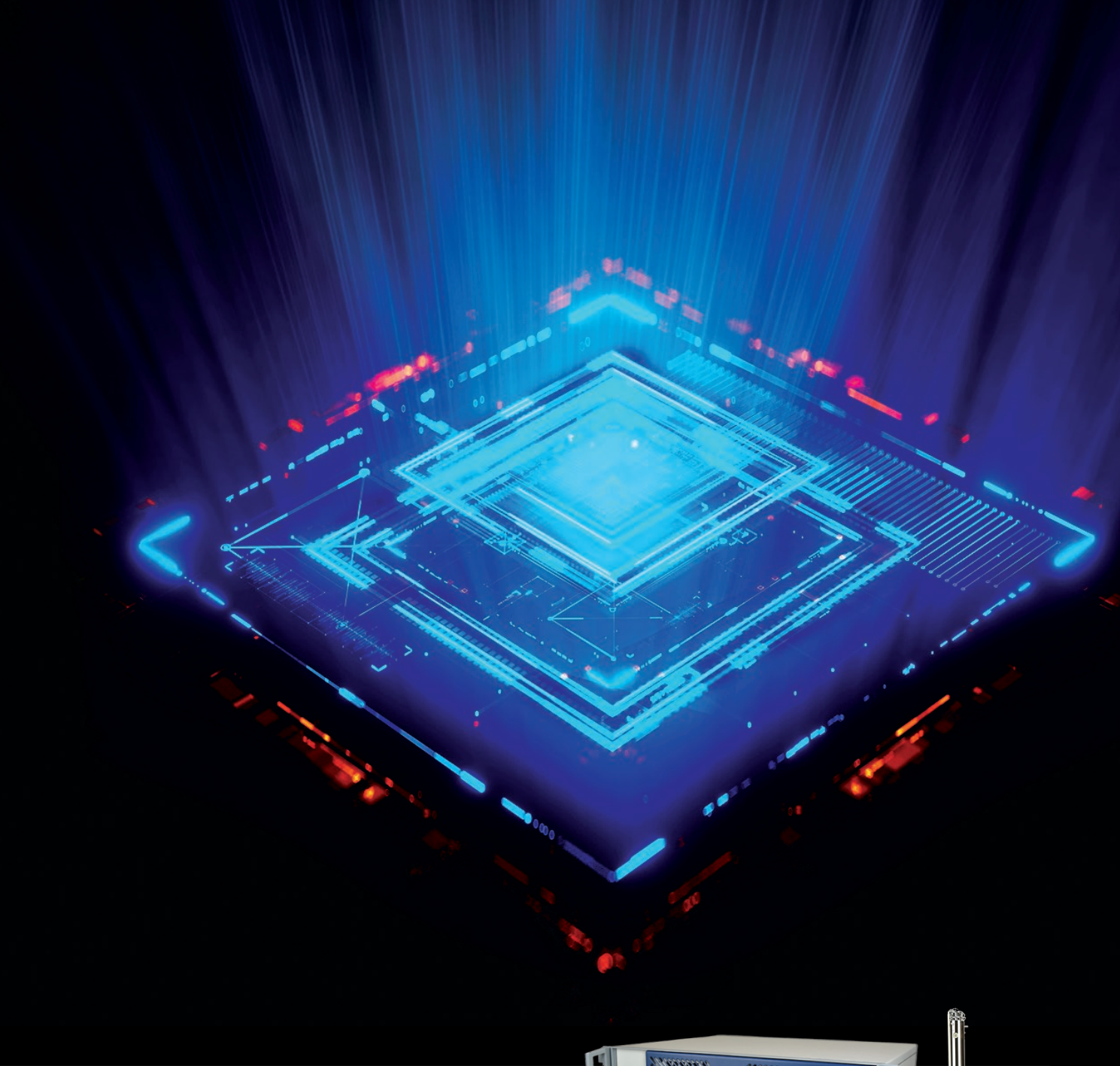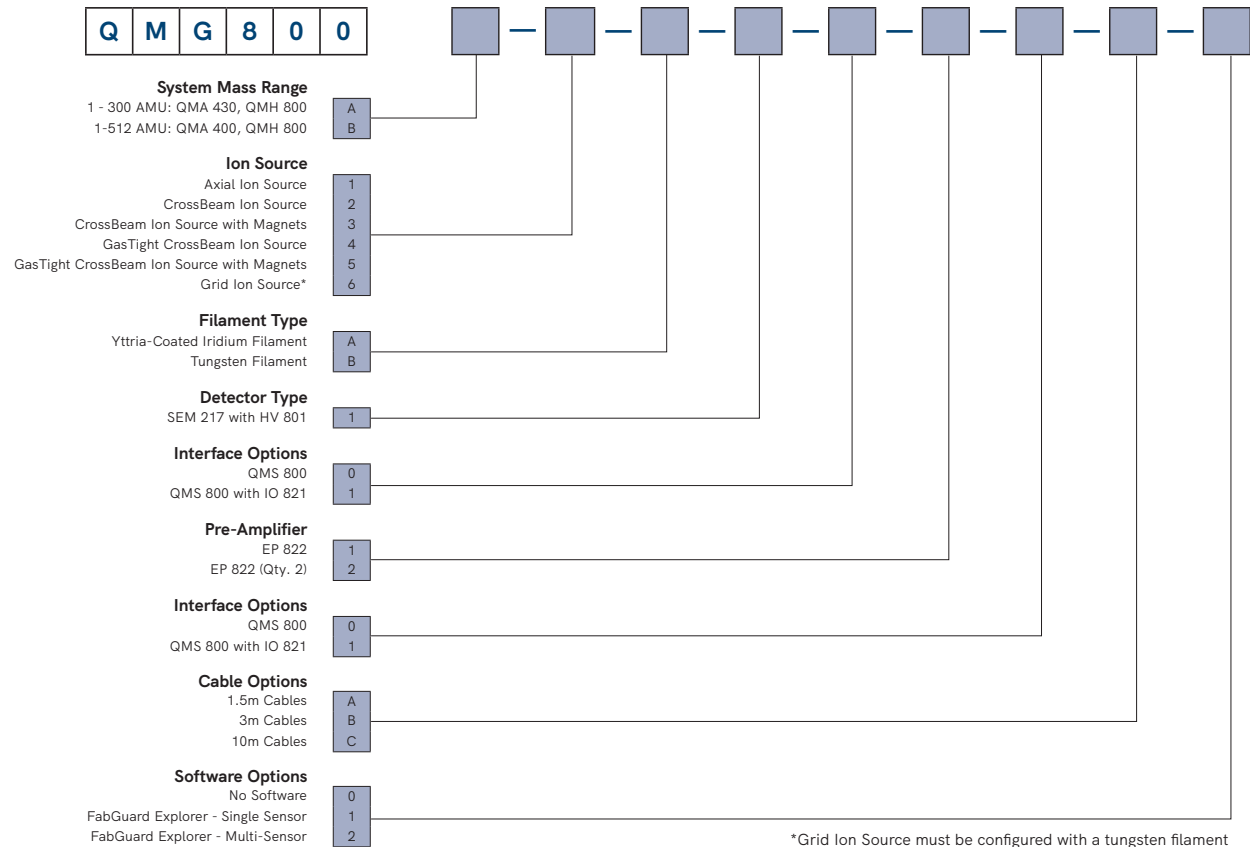

QMG 800 Analytical Mass Spectrometer

Unmatched Sensitivity, Precision, and Analytical Capabilities



Elevate Your Process to New Heights of Innovation

The INFICON QMG 800 is the next-generation analytical mass spectrometer, succeeding QMG 700. QMG 800 remains the market-leading Residual Gas Analyzer (RGA) for analytical process monitoring.


Recognizing customers' unique needs, INFICON built the QMG 800 to allow for more flexibility to meet specific application requirements while maintaining industry-leading measurement speed and sensitivity. QMG 800 is designed to offer enhanced performance and reliability for a wide range of analytical mass spectrometry applications. It is the ideal choice for precise and accurate analysis in contamination monitoring processes, ensuring the highest levels of purity and quality control in various industries. Users will appreciate the digital advances of QMG 800, which incorporates cutting-edge technologies to enhance performance and experience.

With more than 50 years of pioneering excellence in analytical mass spectrometry, INFICON sets a new standard for your gas analysis investments with the QMG 800. Experience unparalleled accuracy and precision, backed by our global expertise and commitment to innovation.

ADVANTAGES AT A GLANCE

- ✓ Advanced Digital Platform
- ✓ Seamless OEM Integration
- ✓ Unparalleled Diagnostics — fully compatible with INFICON FabGuard® software
- ✓ 90° Off-Axis SEM — to minimize noise from radicals and non-ionized particles
- ✓ Single Ion Detection
- ✓ Industry-Leading Measurement

Five ion source options

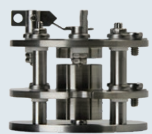

AXIAL

Open design with a single filament; higher sensitivity and good linearity

Dimensions: 26mm H

CROSS BEAM

Open design with dual filament; reacts quickly to change in gas composition and has a longer lifetime



Dimensions: 35.5mm H



CROSS BEAM GAS TIGHT

Closed design with dual filament; in addition to standard benefits, has a sealed ionization area

Dimensions: 48mm H

CROSS BEAM WITH MAGNETS

Allow improved focus of ions to improve sensitivity



Dimensions: 35.5mm H



CROSS BEAM GAS TIGHT WITH MAGNETS

Allow improved focus of ions to improve sensitivity

Dimensions: 48mm H


Sensor | QMA 4x0

QUADRUPOLE

- Allows scanning across the system's mass range.
- 300 AMU or 512 AMU mass range option, utilizing Stainless Steel and Molybdenum, respectively

TYPICAL APPLICATIONS

- Contamination Detection
- Gas Purity Monitor
- End-Point Detection
- Molecular Beam Epitaxy Monitor

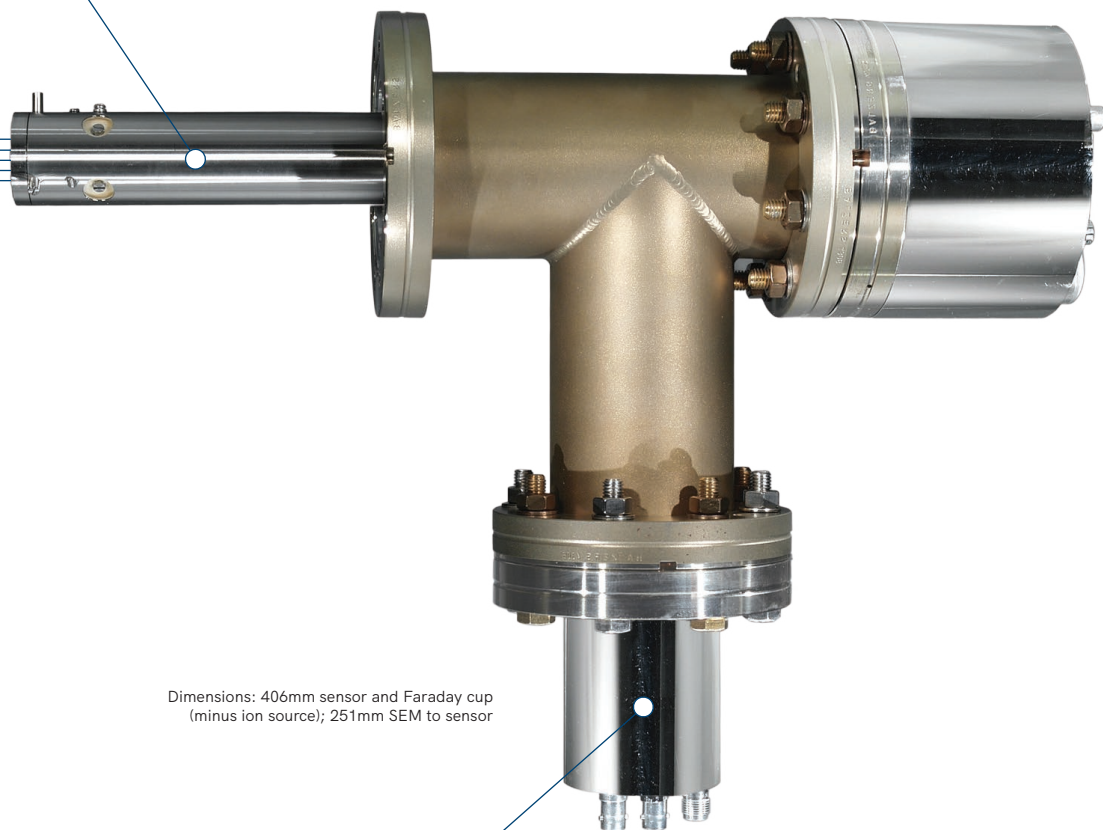

Dimensions: 406mm sensor and Faraday cup (minus ion source); 251mm SEM to sensor

SEM 217

- 90-degree off-axis electron multiplier
- 17 channels to minimize noise



EP 822

- Electrometer preamplifier
- Takes signals from the sensor and amplifies them for readout by the system controller



RF Generator | QMH 800

CUSTOM TRANSFORMER

Designed and built for this product, the customer transformer is the heart of the RF generator, allowing all mechanical and electrical functionality.

PERFORMANCE SPEC

The latest generation boasts improved temperature stability with a less than 0.01 AMU drift per degree Celsius.



Dimensions: 302 mm L x 193 mm W x 229 mm H

INTERNAL OVEN

The latest generation design improves stability across the mass spectrum.

DIGITAL PLATFORM

The new digital platform designed into the RF generator allows for a digital tune functionality.

FABGUARD — DATA COLLECTION AND ANALYSIS

When integrated with the INFICON FabGuard software suite, QMG 800 becomes a powerful process monitoring and diagnostics tool which can be used for:

- Advanced process control (endpoint detection)
- Run-by-run and real-time fault detection and classification
- Statistical process control (SPC)
- Maintenance and decision support with FabRecover

